


Famed LM-01.4



COMFORTABLE, ROBUST AND FUNCTIONAL MATERNITY BED WITH ELECTRICAL CONTROL

Maternity bed Famed LM-01.4 is designed for obstetrics and gynaecology wards patients, both within pre-labour, labour and puerperium period. Maternity bed can be easily transformed into a delivery chair ensuring comfortable position during labour for a future mother. The compact design favors the use of this bed in smaller rooms. Famed LM-01.4 offers continuous, reliable operation thanks to its built-in battery, which is perfect for emergency situations (several cycles).

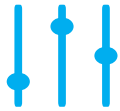

CONTROL OPTIONS

The bed can be controlled using panels in side handrails, which enable control of height, back-rest, allow for turning the lights on and off and alarm activation. Additional functionality is provided by the remote control, which in addition to height and angle adjustment back and thigh segments has also a pre-programmed position for getting down from the bed and examination position and CPR. Easy control is also ensured by foot control.



EXTENSIVE ACCESSORIES

Extensive accessories allow for adopting many labour positions and facilitate performing specialist labour-related, proctoscopy, gynaecology and palpatory examinations. Basic equipment includes arm holders, feet support, knees support and stainless steel bowl. Accessories can be installed on universal rails. The bed is equipped with a set of mattresses secured against penetration of fluids.



SAFETY

Stable handrails are mounted to the backrest and can be lowered if necessary. Polyurethane pads on the hand grips prevent slipping of the hands and provide secure support during labour. All plastic and varnished elements have been made with the use of antibacterial agents in the pSilver™ technology.



Length 2170 mm | Total width 1080 mm | Minimum height 615 mm | Maximum height 985 mm | Back rest inclination angle 70° | Thighs segment adjustment 17° | Trendelenburg position 12° | Reverse-Trendelenburg position 5° | Maximum load 230 kg





Antibacterial, phtalate-free upholstery



Yellow
1F



Light blue
2F



Blue
3F



Turquoise
4F



Light grey
5F



Grey
6F



Red
8F



Black
7F