



# K2k

## MANUAL BED

Registration code: SK-C1



**SAIKANG** / Medical




**Medspares  
Pacific**  
EQUIPMENT  
SPARES  
ACCESSORIES

Medspares Pacific Ltd  
PO Box 2070 Raumati Beach 5032  
Wellington - New Zealand  
P: +64 4 299 7610 - F: +64 4 299 2784  
e: sales@medspares.co.nz - www.medspares.co.nz

Healthier  
your life

# K2k

## FEATURES AND HIGHLIGHTS

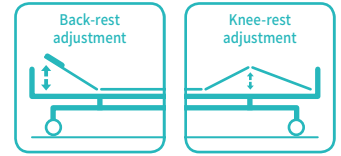


### Technical parameters

▪ External size (LxW):	2115x980 mm
▪ Mattress platform:	1925x820 mm
▪ Safe working load:	200 kg
▪ Height adjustment:	490mm
▪ Back-rest adjustment:	0-75° ( $\pm 10^\circ$ )
▪ Knee-rest adjustment:	0-40° ( $\pm 10^\circ$ )

### Technical configuration

<input checked="" type="checkbox"/> Manual Cranks	2pcs
<input checked="" type="checkbox"/> 5" Covered Castors	4pcs
<input checked="" type="checkbox"/> Foldable Side Rail	1set
<input checked="" type="checkbox"/> ABS Bed Ends with Safe Lock	1set
<input checked="" type="checkbox"/> IV Pole	1set
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs
<input checked="" type="checkbox"/> Mattress Holder	2pcs

### Optional color (MOQ100pcs- color changing)

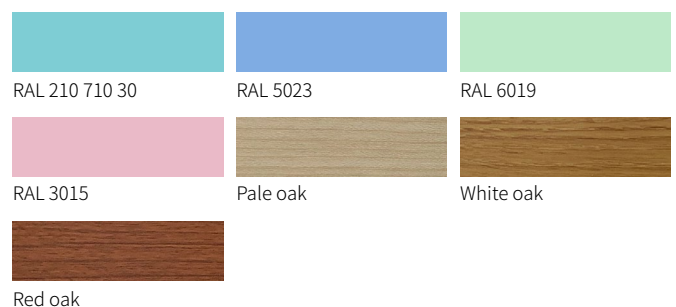

# QUALITY CONTROL

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating, efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



## Design

The design is both modern and practical. The head and foot boards feature a beautiful and elegant shape. The use of a hollow protective wrap angle, combined with a barrel-type bumper, makes the product instantly recognizable.



## Steels | welding

Baosteel from Fortune Global 500. Panasonic Robotic ensure 360° full smooth welding.



## Plastic | Painting

Environmental plastic can be sterilized at a temperature of 100°C and has a tensile strength of up to 30MPa.

11 Process epoxy painting with a thickness of 0.12mm and a brightness of 60°, which is capable of resisting impacts of up to 50kg.

## Accessories



### Head & foot board

The head and foot boards feature "integrated protective corners" with hollow design, which are combined with barrel-type bumper wheels that have a diameter of 60mm.

This design not only minimizes the overall space occupied by the bed but also provides protection.

Additionally, the head and foot boards are effortless to disassemble and clean.




### Castors

5" covered casters.

PR material, walks 30KM with no wear and tear. Precision-designed anti-winding shell with no screws.

Pass Dynamic test: withstand a weight of 120kg, walk 30KM, and pass through 500 obstacles without any deformation or damage.

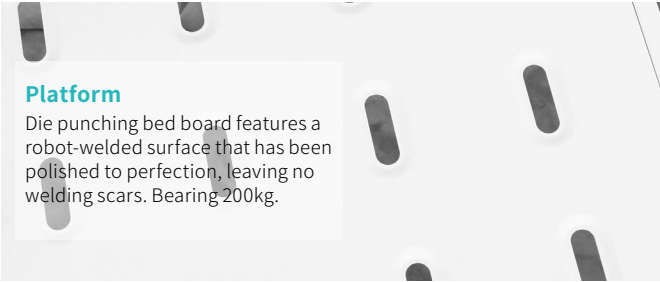


**Platform**

Die punching bed board features a robot-welded surface that has been polished to perfection, leaving no welding scars. Bearing 200kg.



**Siderail**

The guardrail is made of six-pole folding aluminum alloy, with stainless steel poles and a plastic handle.

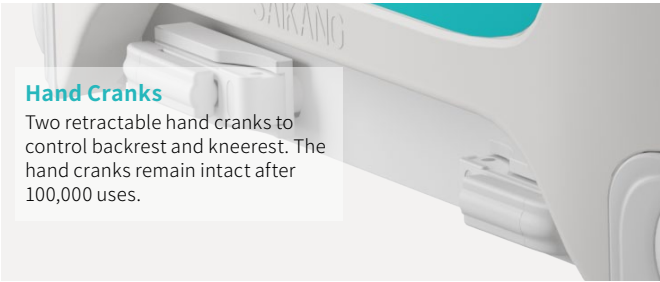


**Hand Cranks**

Two retractable hand cranks to control backrest and kneerest. The hand cranks remain intact after 100,000 uses.



**Mattress Retainer**

The dual-function mattress retainer serves to not only prevent scratches on the bed frame's paint but also retracts the mattress to its place, and prevents mattress sliding.

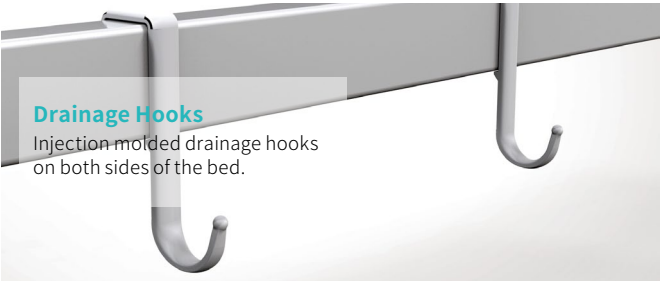

**Drainage Hooks**

Injection molded drainage hooks on both sides of the bed.



**IV Pole**

The infusion pole is made of stainless steel and has four plastic hooks that can bear up to 15 kg of weight.



Safe working load

200 kg



Back-rest adjustment

75°



Knee-rest adjustment

40°